

# FINAL REPORT

# A Pharmacokinetic Study of OCT-598 Following Single Intravenous and Oral Gavage Administration to Sprague-Dawley Rats

**Study Number: 0040523195** 

Nonclinical Research Center QuBEST BIO Co., Ltd.

# **COMPLIANCE STATEMENT AND APPROVAL**

I, the undersigned, hereby declare that the work was performed under my supervision and that the report represents a true and accurate record of the results obtained.

This study was performed in accordance with the agreed protocol and with Standard Operating Procedures, unless otherwise stated, and the study objectives were achieved.

This study is not within the scope of regulations governing the conduct of nonclinical laboratory studies and is not intended to comply with such regulations.

| Study Director: | Nov-24-2023                                  |           |
|-----------------|--|-----------|
|                 | Cheolmin Park, BS                            | Date      |
|                 | Nonclinical Research Center, QuBEST BIO C    | Co., Ltd. |
|                 | TEL: +82-31-5181-8743, FAX: +82-31-5181-8701 |           |
|                 | E-mail: pcm@qubest.co.kr                     |           |

Management:

Soohyeon Kim, MS

QuBEST BIO Co., Ltd.

TEL: +82-31-5181-8777, FAX: +82-31-706-2996

E-mail: shkim@qubest.co.kr

# TABLE OF CONTENTS

| CO  | MPLI/  | ANCE STATEMENT AND APPROVAL   | 2   |
|-----|--------|---|-----|
| 1.  | UNIT   | S AND ABBREVIATIONS   | 5   |
| 2.  | SUM    | MARY  | 6   |
| 3.  | STUI   | DY INFORMATION  | 8   |
|     | 3.1    | Objectives  | 8   |
|     | 3.2    | Regulatory Test Guidelines  | 8   |
|     | 3.3    | Sponsor   | 8   |
|     | 3.4    | Test Facility   | 8   |
|     | 3.5    | Study Schedule  | 8   |
|     | 3.6    | Responsible Personnel   | 9   |
| 4.  | MAT    | ERIALS AND METHODS  | 9   |
|     | 4.1    | Test Article and Vehicle Information  | 9   |
|     | 4.2    | Preparation of Dose Formulations  | 10  |
|     | 4.3    | Test System   | 10  |
|     | 4.4    | Study Method  | 12  |
|     | 4.5    | In-life Observations and Examinations   | 13  |
|     | 4.6    | Sample Collection   | 13  |
|     | 4.7    | Bioanalysis   | 13  |
|     | 4.8    | Data Processing and Pharmacokinetic Analysis                                    | 14  |
| 5.  | RETE   | ENTION OF RECORDS, SAMPLES AND SPECIMENS  | 14  |
| 6.  | RESU   | JLTS AND DISCUSSION   | 14  |
| 7.  | CON    | CLUSIONS  | 16  |
| 8.  | REFE   | ERENCES   | 17  |
| 9.  | FIGU   | RES   | 18  |
|     | Figure | e 1. Mean (±SD) Plasma Concentrations-Time Profiles of OCT-598 Following Single |     |
|     | Intrav | enous Administration to Male Rats at 3 mg/kg                                    | 18  |
|     | Figure | e 2. Mean (±SD) Plasma Concentrations-Time Profiles of OCT-598 Following Single |     |
|     | Oral ( | Gavage Administration to Fasted Male Rats at 3, 10 and 30 mg/kg                 | 19  |
|     | Figure | e 3. Mean (±SD) Plasma Concentrations-Time Profiles of OCT-598 Following Single |     |
|     | Oral ( | Gavage Administration to Fasted Male and Female Rats at 10 mg/kg                | 20  |
|     | Figure | e 4. Mean (±SD) Plasma Concentrations-Time Profiles of OCT-598 Following Single |     |
|     | Oral O | Gavage Administration to Fasted and Fed Male Rats at 10 mg/kg                   | 21  |
| 10. | TABI   | JE  | 22  |
|     | Table  | 1. Clinical Observations and Body Weights (Group Summary)                       | .22 |

| 11. | APPENDICES   | . 23 |
|-----|--|------|
|     | Appendix 1. Individual Clinical Observations and Body Weights                          | . 23 |
|     | Appendix 2. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 |      |
|     | Following Single Intravenous Administration to Fed Male Rats at 3 mg/kg                | . 25 |
|     | Appendix 3. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 |      |
|     | Following Single Oral Gavage Administration to Fasted Male Rats at 3 mg/kg             | 26   |
|     | Appendix 4. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 |      |
|     | Following Single Oral Gavage Administration to Fasted Female Rats at 10 mg/kg          | 27   |
|     | Appendix 5. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 |      |
|     | Following Single Oral Gavage Administration to Fasted Male Rats at 10 mg/kg            | 28   |
|     | Appendix 6. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 |      |
|     | Following Single Oral Gavage Administration to Fasted Male Rats at 30 mg/kg            | 29   |
|     | Appendix 7. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 |      |
|     | Following Single Oral Gavage Administration to Fed Male Rats at 10 mg/kg               | 30   |
|     | Appendix 8. Certificate of Analysis (COA)  | .31  |
|     | Appendix 9. Analytical Method  | 33   |

# 1. UNITS AND ABBREVIATIONS

The following lists of codes, abbreviations and units are used by QuBEST BIO Co., Ltd. All of the abbreviations listed on these pages may not be applicable to this Protocol.

# GLOSSARY OF ABBREVIATIONS AND DEFINITION OF TERMS

| %       | Percent                              | BLQ      | Below the limit of quantification |
|---------|--------------------------------------|----------|-----------------------------------|
| °C      | Degree Celsius                       | Conc.    | Concentration                     |
| kg      | Kilogram                             | LC-MS/MS | HPLC-tandem mass spectrometry     |
| g       | Gram                                 | LLOQ     | Lower limit of quantification     |
| mg      | Milligram (10 <sup>-3</sup> grams)   | NA       | Not applicable                    |
| μg      | Microgram (10 <sup>-6</sup> grams)   | ND       | Not detected/Not determined       |
| ng      | Nanogram (10 <sup>-9</sup> grams)    | No.      | Number                            |
| L       | Liter                                | SD       | Standard deviation                |
| mL      | Milliliter (10 <sup>-3</sup> liters) | TBD      | To be determined                  |
| $\mu L$ | Microliter (10 <sup>-6</sup> liters) | ULOQ     | Upper limit of quantification     |
| min     | Minute                               | CV (%)   | Percent coefficient of variation  |
| hr      | Hour                                 |          |                                   |

# PHARMACOKINETIC ABBREVIATIONS

| AUClast      | Area under the curve from the time of dosing to the last measurable concentration   |  |  |
|--------------|---|--|--|
| AUCinf       | Area under the curve from dosing time extrapolated to infinity  |  |  |
| C0           | Initial concentration   |  |  |
| Cmax         | The maximum measured concentration  |  |  |
| CL           | Total clearance after intravenous administration  |  |  |
| CL/F         | Apparent total clearance after extravascular administration   |  |  |
| Tmax         | The time to reach maximum concentration   |  |  |
| t1/2         | Terminal half-life  |  |  |
| Rsq_adjusted | Goodness of fit statistic for the terminal elimination phase, adjusted for the number of points used in the estimation of Lambda Z.                       |  |  |
| Vdss         | Volume of distribution at steady-state after intravenous administration   |  |  |
| Vd/F         | Apparent volume of distribution after extravascular administration  |  |  |
| Vz           | Volume of distribution based on the terminal phase  |  |  |
| %AUCexp      | Percentage of the AUCinf that is contributed by the extrapolation from the last sampling time to infinity   |  |  |
| BA           | Bioavailability (%) = $\frac{\text{AUC}extravescular}{\text{AUC}intravenous} \times \frac{\text{Dose}extravescular}{\text{Dose}extravescular} \times 100$ |  |  |

#### 2. SUMMARY

The objective of this study was to evaluate the pharmacokinetic profiles of OCT-598 following single intravenous and oral gavage administration to Sprague-Dawley rats.

OCT-598 was administered to male and female rats in fasted or fed conditions at 3, 10 and 30 mg/kg via oral gavage and 3 mg/kg via intravenous route. The blood samples were collected at the scheduled times up to 24 hours post-dose. Determination of OCT-598 in plasma was analyzed by UHPLC-MS/MS. Pharmacokinetic parameters were calculated by noncompartmental analysis using Phoenix® WinNonlin® (Ver. 8.4, Certara).

Following single intravenous administration of OCT-598 to male rats at 3 mg/kg, plasma concentration of OCT-598 decreased with a terminal half-life (t1/2) of 4.08 hours. Total clearance (CL) of OCT-598 was 1099.82 mL/hr/kg, slower than hepatic blood flow in rats. Volume of distribution at steady-state (Vdss) of OCT-598 was 2466.54 mL/kg, larger than total body water indicating that OCT-598 was more distributed into tissue than circulating blood.

Following single oral gavage administration of OCT-598 to fasted male rats at 3, 10 and 30 mg/kg, OCT-598 was absorbed with a median Tmax of 0.5 hours. After reaching Cmax, plasma concentration of OCT-598 decreased with a terminal half-life (t1/2) ranged from 2.64 to 3.49 hours. Cmax ranged from 671.67 to 6411.18 ng/mL and AUClast ranged from 1692.91 to 25456.48 ng·hr/mL. Systemic exposure (Cmax and AUClast) of OCT-598 increased more than dose-proportionally. Specifically, as the dose increased in a ratio of 1.0 : 3.3 : 10.0, the Cmax increased in a ratio of 1.0 : 5.1 : 9.5 and the AUClast increased in a ratio of 1.0 : 4.7 : 15.0.

The oral bioavailability of OCT-598 at 3 mg/kg was 62.41%.

Following single oral gavage administration of OCT-598 to fasted female rats at 10 mg/kg, OCT-598 was absorbed with a median Tmax of 0.5 hours. After reaching Cmax, plasma concentration of OCT-598 decreased with a terminal half-life (t1/2) of 3.81 hours. Cmax and AUClast were 3701.88 ng/mL and 6106.66 ng·hr/mL, respectively.

Following single oral gavage administration of OCT-598 to fed male rats at 10 mg/kg, OCT-598 was absorbed with a median Tmax of 1 hour. After reaching Cmax, plasma concentration of OCT-598 decreased with a terminal half-life (t1/2) of 4.99 hours. Cmax and AUClast were 976.03 ng/mL and 6069.72 ng·hr/mL, respectively.

No notable gender difference (<2-fold) was observed.

The food effect was observed in terms of Tmax and Cmax. Tmax was delayed by approximately 0.5 hour by food. Compared with the fasting group, Cmax was altered significantly with a decrease of approximately 72% and AUClast did not alter significantly (<2-fold).

#### 3. STUDY INFORMATION

# 3.1 Objectives

The objective of this study was to evaluate the pharmacokinetic profiles of OCT-598 following single intravenous (bolus) and oral gavage administration to Sprague-Dawley rats.

#### 3.2 Regulatory Test Guidelines

This is non-GLP study and was not conducted in compliance with GLP regulation; however, it was conducted according to the protocol and protocol amendment (if any) approved by the Sponsor and SOPs of Nonclinical Research Center, QuBEST BIO.

## 3.3 Sponsor

Oscotec Inc.

9th Floor, Building A, 700 Daewangpangyo-ro, Bundang-gu, Seongnam-si, Gyeonggi-do, 13488, Republic of Korea

TEL: +82-31-628-7631, FAX: +82-31-628-7668

# 3.4 Test Facility

Nonclinical Research Center, QuBEST BIO Co., Ltd.

#301, Daewoo Frontier Valley I, 16-25, Dongbaekjungang-ro 16beon-gil, Giheung-gu, Yongin-si, Gyeonggi-do, 17715, Republic of Korea

TEL: +82-31-5181-8700, FAX: +82-31-5181-8701

#### 3.5 Study Schedule

| Study Initiation:           | September 14, 2023                      |  |
|-----------------------------|---|--|
| Experimental Start:         | September 14, 2023                      |  |
| Administration:             | September 20, 2023                      |  |
| Sample Collection:          | September 20, 2023 - September 21, 2023 |  |
| Pharmacokinetic Analysis:   | September 21, 2023 - September 27, 2023 |  |
| Submission of Final Report: | November 24, 2023                       |  |

#### 3.6 Responsible Personnel

Storage/Preparation of the Test Article:

Animal Experiment:

Geumjin Choi, BS, KLAT

Jiyu Oh, BS, KLAT

Junhui Lim, BS

Daeyeon In, BS

Bioanalysis:

Sungyong Choi, MS

Jaeyoung Seo, BS

Internal Scientific Review:

Soohyeon Kim, MS

Archives:

Myeongseok Gu, BS

# 4. MATERIALS AND METHODS

# 4.1 Test Article and Vehicle Information

#### 4.1.1 Test Article

| CPo131731-01-08-01   |
|--|
| White solid  |
| 97.3%  |
| 507.62   |
| 507.19   |
| Not applied  |
| May 26, 2024   |
| Refrigerate (2-8°C), protected from light and dehumidification |
| Oscotec Inc.   |
| Returned to the Sponsor  |
|  |

#### 4.1.2 Vehicle

| Name:               | 10% SBE-β-CD in 50 mM Potassium phosphate buffer (pH 8.0) |
|---------------------|---|
| Storage Conditions: | Room temperature  |
| Manufacturer:       | Nonclinical Research Center, QuBEST BIO                   |

## 4.1.2.1 Vehicle Component

| Name   | Batch/Lot Number | Supplier              |
|--|------------------|-----------------------|
| SBE-β-CD (Sulfobutylether-beta-cyclodextrin) | OC15791801       | BIOSYNTH (Carbosynth) |
| Potassium phosphate dibasic                  | SLCK8159         | Sigma-Aldrich         |
| Potassium phosphate monobasic                | NKCR5040         | Sigma-Aldrich         |
| DW (Distilled water)                         | S5X8B21          | DAI HAN PHARM         |

#### 4.2 Preparation of Dose Formulations

The dose formulations were freshly prepared in clean bench on the day of each dosing according to the procedure provided by the Sponsor and a correction factor was not applied to prepare the dose formulations.

- 1) Required amount of test article was weighed and transferred to appropriate container.
- 2) 10% SBE-β-CD in 50 mM potassium phosphate buffer (approximately 80% of total volume) was added into the container and stirred/sonicated until clear solution by visual check (The dosing formulations were sonicated at least 30 minutes).
- 3) Finally, 10% SBE-β-CD in 50 mM potassium phosphate buffer (approximately 20% of total volume) was added into the container to the final formulation volume and stirred/vortexed until clear solution by visual check.

During the dosing, formulations were handled at room temperature, and the remaining formulations were discarded.

#### 4.2.1 Analysis of Dose Formulations

Analysis of homogeneity, concentration verification and stability of the dose formulations were not performed in the Test Facility.

#### 4.3 Test System

#### 4.3.1 Animal Information

| Species and Strain:                       | Specific pathogen free (SPF) rat, Sprague-Dawley [NTacSam:SD]     |  |
|---|---|--|
| Breeder/Supplier:                         | SAMTAKOBIOKOREA, Inc.   |  |
|   | (105, Seorang-ro, Osan-si, Gyeonggi-do, 18100, Republic of Korea) |  |
| Sex and Number of Animals Ordered:        | Males, 33 and Females, 7  |  |
| Age at Receipt:                           | Approximately 7 weeks   |  |
| Age at the Initiation of Dosing:          | Approximately 8 weeks   |  |
| Weight Range at the Initiation of Dosing: | 159.9 - 294.5 g   |  |

#### 4.3.2 Justification for Selection

The Sprague-Dawley rats were chosen in this study because they were widely used in the pharmacokinetic study of drugs. The number of animals to be used in this study was considered the minimum number required to evaluate the pharmacokinetics of the test article.

#### 4.3.3 Identification

All animals were individually identified by tail marking method using indelible red marker during acclimation period. During the study, each animal was identified by tail marking method using indelible black marker and a cage label card displaying the study number, group, and animal number.

#### 4.3.4 Animal Welfare

The Protocol and procedures involving the care and use of animals in this study was reviewed and approved by QuBEST BIO IACUC prior to conduct (Approval No.: QBIACUC-A23195). During the study, the care and use of animals were conducted in accordance with all applicable guidelines of Animal Protection Law.

## 4.3.5 Husbandry and Environmental Conditions

This study was performed in the barrier animal facility area of Nonclinical Research Center, QuBEST BIO.

| HVAC Conditions:           | 100% HEPA-filtered air, at least 10~20 air changes/hr           |  |
|----------------------------|---|--|
| Temperature and Humidity:  | 20.7 - 21.7°C, 57.3 - 61.3% (relative humidity)                 |  |
| Light Cycle and Intensity: | 12 hours light and 12 hours dark (on 08:00-20:00, except during |  |
|                            | designated procedures in the Protocol), intensity 150-300 Lux   |  |

The animals were housed in polycarbonate cages [260W x 420L x 180H (mm)] with irradiated bedding material, up to three per cage. Animal room and each cage/bottle were cleaned at regular intervals per SOPs of Nonclinical Research Center, QuBEST BIO.

#### 4.3.6 Food and Water

All animals were offered to standard irradiated pelleted commercial laboratory diet (Purina Rodent Chow 38057, Republic of Korea) and water during the study. Nutritional components and environmental contaminants in the food were analyzed routinely by the manufacturer. No contaminants were reasonably expected to be present that would interfere with the objectives of the study; therefore, no testing was conducted as part of the study. According to SOPs of Nonclinical Research Center in QuBEST BIO, the drinking water was analyzed for contaminants semiannually by an independent laboratory. Results of the analyses of food and water are on file at the Test Facility.

#### 4.3.7 Randomization

On the last day of acclimation period, all animals were weighed, evaluated for general health and suitability of testing and those considered suitable for the study was released to the study. The animals were randomly assigned to study groups by a stratified randomization scheme designed to achieve similar group mean body weights. Prior to the initiation of dosing and/or immediately after dosing, any

assigned animals considered unsuitable for use and/or accidental events were replaced from the remaining unassigned animals.

#### 4.3.8 Food Restriction

| Dosing Route (Group)   | Fasting Time           | Food Supply                     |
|------------------------|------------------------|---------------------------------|
| Intravenous (G1)       | Not applicable         | Ad libitum                      |
| Oral gavage (G2 to G5) | Approximately 16 hours | Approximately 4 hours post-dose |
| Oral gavage (G6)       | Not applicable         | Ad libitum                      |

#### 4.4 Study Method

#### 4.4.1 Group Assignment

| Canada        | Tuestment | Dosing<br>Route | Dose<br>Level | Dose<br>Volume<br>(mL/kg) | No. of  | Animal ID |           | Food        |
|---------------|-----------|-----------------|---------------|---------------------------|---------|-----------|-----------|-------------|
| Group Treatme | Treatment |                 | (mg/kg)       |                           | Animals | Males     | Females   | Restriction |
| G1            |           | IV              | 3             | 1                         | 3       | M1 - M6   |           | Fed         |
| G2            | -         |                 | 3             |                           | 3       | M7 - M12  |           |             |
| G3            | OCT 500   | OCT 509         | 10            |                           | 3       |           | F13 - F18 | Footod      |
| G4            | OCT-598   | PO              | 10            | 5                         | 3       | M19 - M24 |           | Fasted      |
| G5            |           |                 | 30            | -                         | 3       | M25 - M30 |           |             |
| G6            | -         |                 | 10            | -                         | 3       | M31 - M36 |           | Fed         |

#### 4.4.2 Justification of Administration Routes

The oral gavage administration route was selected because this is the intended route of human exposure. The intravenous administration was selected to investigate the total clearance, volume of distribution and oral bioavailability of test article.

#### 4.4.3 Justification of Dose Levels

Dose levels were selected by the Sponsor.

#### 4.4.4 Justification of Dose Volumes

Dose volumes were selected by the Sponsor.

#### 4.4.5 Method of Administration

#### 4.4.5.1 Intravenous Administration

The dose formulation was administered intravenously into tail vein via a disposable syringe with 26-gauge needle. Individual dose volumes were based on the most recent body weights.

#### 4.4.5.2 Oral Gavage Administration

The dose formulation was administered orally via a syringe equipped with a disposable feeding tube [straight, 15-gauge (3.0 mm bulb diameter x 100 mm length)]. Individual dose volumes were based on the most recent body weights.

#### 4.4.6 Dosing Regimen

The dose formulation was administered once. The day of dosing was designated as Day 1.

## 4.5 In-life Observations and Examinations

#### 4.5.1 Mortality and Clinical Observations

All animals were observed twice daily for mortality, abnormalities, and signs of pain and stress, once in the morning and once in the afternoon.

#### 4.5.2 Body Weights

Individual body weights were measured prior to dosing for all animals.

#### 4.6 Sample Collection

#### 4.6.1 Blood Collection

| Sampling Method:   | Composite sampling (3 rats/time point)   |
|--------------------|--|
| Sampling Schedule: | IV: 0.083, 0.25, 0.5, 1, 2, 4, 8 and 24 hours post-dose (total 8 time points) PO: 0.5, 1, 2, 4, 8, 12 and 24 hours post-dose (total 7 time points) |
| Sample Volume:     | Approximately 1 mL/time point  |
| Sampling Site:     | Jugular vein   |
| Anticoagulant:     | Sodium heparin   |
| Sample Handling:   | Keep samples chilled (wet ice, as appropriate) during collection and processing  |

#### 4.6.2 Plasma Sample Processing

| Centrifugation:    | 12,000 rpm for approximately 2 minutes at 4°C                        |
|--------------------|--|
| Aliquots:          | The maximum amount of plasma was recovered and placed in uniquely    |
|                    | labeled micro tubes. Plasma was divided into two approximately equal |
|                    | aliquots.  |
| Label Information: | Study group, animal number and time point                            |

#### 4.6.3 Sample Storage

Plasma samples were stored frozen (below -70°C) unless otherwise specified in a validated method.

#### 4.7 Bioanalysis

The concentrations of OCT-598 in plasma was analyzed at Bioanalysis Center, QuBEST BIO using a qualified UHPLC-MS/MS method.

#### 4.8 Data Processing and Pharmacokinetic Analysis

Pharmacokinetic parameters were calculated at Nonclinical Research Center, QuBEST BIO. Non-compartmental analysis of OCT-598 concentrations in plasma was performed by using the Phoenix® WinNonlin® software (Ver. 8.4, Certara).

Plasma concentrations below the LLOQ were replaced with zero for calculation purposes.

The area under the concentration versus time curve (AUC) was calculated using linear trapezoidal method with linear interpolation. The slope of the terminal elimination phase was determined using the best fit for  $\lambda z$  and log-linear regression on the unweighted concentration data.

#### 5. RETENTION OF RECORDS, SAMPLES AND SPECIMENS

The study-specific raw data, documents, correspondences, Protocol and Final Report will be archived at the test facility for six months after issuance of Final Report. Remained test or reference articles, specimens and study samples will be retained under appropriate condition by the issuance of Final Report but not later than six months after submission of Draft Report. The Sponsor will be informed the list of the specimens and study samples remained prior to discard.

#### 6. RESULTS AND DISCUSSION

During the in-life period, no test article-related abnormality was observed in all study animals (Table 1 and Appendix 1). The mean pharmacokinetic parameters and concentration vs. time profiles of OCT-598 are presented in Text Table 1 to Text Table 7 and Figure 1 to Figure 4. Individual pharmacokinetic parameters of OCT-598 are presented in Appendix 2 to Appendix 7.

Text Table 1. Mean Pharmacokinetic Parameters of OCT-598 Following Single Intravenous Administration to Male Rats at 3 mg/kg

| DV Damanatana      | Male (Fed) |  |
|--------------------|------------|--|
| PK Parameters      | 3 mg/kg    |  |
| C0 (ng/mL)         | 3103.60    |  |
| AUClast (ng·hr/mL) | 2712.40    |  |
| AUCinf (ng·hr/mL)  | 2727.72    |  |
| t1/2 (hr)          | 4.08       |  |
| CL (mL/hr/kg)      | 1099.82    |  |
| Vdss (mL/kg)       | 2466.54    |  |
| Vz (mL/kg)         | 6479.99    |  |
| Rsq_adjusted       | 0.99       |  |
| %AUCexp (%)        | 0.56       |  |

Text Table 2. Mean Pharmacokinetic Parameters of OCT-598 Following Single Oral Gavage Administration to Fasted Male Rats at 3, 10 and 30 mg/kg

| DIZ D               |         | Male (Fasted) |          |
|---------------------|---------|---------------|----------|
| PK Parameters —     | 3 mg/kg | 10 mg/kg      | 30 mg/kg |
| Cmax (ng/mL)        | 671.67  | 3444.57       | 6411.18  |
| Tmax (hr)           | 0.5     | 0.5           | 0.5      |
| AUClast (ng·hr/mL)  | 1692.91 | 8010.87       | 25456.48 |
| AUCinf (ng·hr/mL)   | 1703.68 | 8027.22       | 25692.25 |
| t1/2 (hr)           | 3.49    | 2.64          | 3.32     |
| CL/F (mL/hr/kg)     | 1760.89 | 1245.76       | 1167.67  |
| Vd/F (mL/kg)        | 8859.37 | 4737.52       | 5591.06  |
| Rsq_adjusted        | 0.98    | 0.98          | 0.96     |
| %AUCexp (%)         | 0.63    | 0.2           | 0.92     |
| Bioavailability (%) | 62.41   | NA            | NA       |

Text Table 3. Mean Pharmacokinetic Parameters of OCT-598 Following Single Oral Gavage Administration to Fasted Male and Female Rats at 10 mg/kg

| PK Parameters —    | Male (Fasted) | Female (Fasted) |
|--------------------|---------------|-----------------|
| PK Parameters      | 10 mg/kg      | 10 mg/kg        |
| Cmax (ng/mL)       | 3444.57       | 3701.88         |
| Tmax (hr)          | 0.5           | 0.5             |
| AUClast (ng·hr/mL) | 8010.87       | 6106.66         |
| AUCinf (ng·hr/mL)  | 8027.22       | 6177.47         |
| t1/2 (hr)          | 2.64          | 3.81            |
| CL/F (mL/hr/kg)    | 1245.76       | 1618.79         |
| Vd/F (mL/kg)       | 4737.52       | 8906.93         |
| Rsq_adjusted       | 0.98          | 0.99            |
| %AUCexp (%)        | 0.2           | 1.15            |

Text Table 4. Mean Pharmacokinetic Parameters of OCT-598 Following Single Oral Gavage Administration to Fasted and Fed Male Rats at 10 mg/kg

| DV Damanastana     | Male (Fasted) | Male (Fed) |
|--------------------|---------------|------------|
| PK Parameters ——   | 10 mg/kg      | 10 mg/kg   |
| Cmax (ng/mL)       | 3444.57       | 976.03     |
| Tmax (hr)          | 0.5           | 1          |
| AUClast (ng·hr/mL) | 8010.87       | 6069.72    |
| AUCinf (ng·hr/mL)  | 8027.22       | 6303.04    |
| t1/2 (hr)          | 2.64          | 4.99       |
| CL/F (mL/hr/kg)    | 1245.76       | 1586.54    |
| Vd/F (mL/kg)       | 4737.52       | 11421.73   |
| Rsq_adjusted       | 0.98          | 0.96       |
| %AUCexp (%)        | 0.2           | 3.7        |

Text Table 5. Dose Proportionality of OCT-598 Following Single Oral Gavage Administration to Fasted Male Rats at 3, 10 and 30 mg/kg

| Ratio         | 3 mg/kg | 10 mg/kg | 30 mg/kg |
|---------------|---------|----------|----------|
| Dose ratio    | 1.0     | 3.3      | 10.0     |
| Cmax ratio    | 1.0     | 5.1      | 9.5      |
| AUClast ratio | 1.0     | 4.7      | 15.0     |

Text Table 6. Gender Difference of OCT-598 Following Single Oral Gavage Administration to Male and Female Rats at 10 mg/kg

| PK Parameters      | Male    | Female  | Ratio <sup>1)</sup> |
|--------------------|---------|---------|---------------------|
| Cmax (ng/mL)       | 3444.57 | 3701.88 | 0.93                |
| AUClast (ng·hr/mL) | 8010.87 | 6106.66 | 1.31                |

<sup>1):</sup> Male/Female

Text Table 7. Food Effect of OCT-598 Following Single Oral Gavage Administration to Male Rats at 10 mg/kg

| PK Parameters      | Fasted  | Fed     | Ratio <sup>1)</sup> |
|--------------------|---------|---------|---------------------|
| Cmax (ng/mL)       | 3444.57 | 976.03  | 0.28                |
| AUClast (ng·hr/mL) | 8010.87 | 6069.72 | 0.76                |

<sup>1):</sup> Fed/Fasted

#### 7. CONCLUSIONS

This study was performed to understand the pharmacokinetic profiles of OCT-598 and to evaluate the dose proportionality, bioavailability, gender difference and food effect by single intravenous or oral gavage administrations. in Sprague-Dawley rats.

Following single intravenous administration of OCT-598 to male rats at 3 mg/kg, total clearance (CL), volume of distribution at steady-state (Vdss) and terminal half-life (t1/2) were 1099.82 mL/hr/kg, 2466.54 mL/kg and 4.08 hours, respectively.

Following single oral gavage administration of OCT-598 to fasted male rats at 3, 10 and 30 mg/kg, systemic exposure (Cmax and AUClast) increased more than dose-proportionally and oral bioavailability of OCT-598 at 3 mg/kg was 62.41%.

No notable gender difference was observed.

The food effect was observed in terms of Tmax and Cmax. Tmax was delayed by approximately 0.5 hour by food. Compared with the fasting group, Cmax was altered significantly with a decrease of approximately 72% and AUClast did not alter significantly (<2-fold).

# 8. REFERENCES

- 1) B. Davies and T. Morris, "Physiological Parameters in Laboratory Animals and Humans", Pharmaceutical Research, Vol. 10, No. 7, 1993, pp. 1093-1095.
- 2) Katya Tsaioun and Steven A. Kates, "ADMET for Medicinal Chemists: A Practical Guide", John Wiley & Sons, 2011, p159.
- 3) Kwon, Younggil, "Handbook of Essential Pharmacokinetics, Pharmacodynamics, and Drug Metabolism for Industrial Scientists", Springer, 2002, p74

# 9. FIGURES

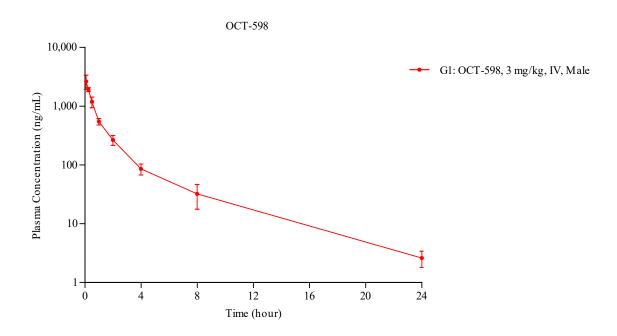


Figure 1. Mean (±SD) Plasma Concentrations-Time Profiles of OCT-598 Following Single Intravenous Administration to Male Rats at 3 mg/kg

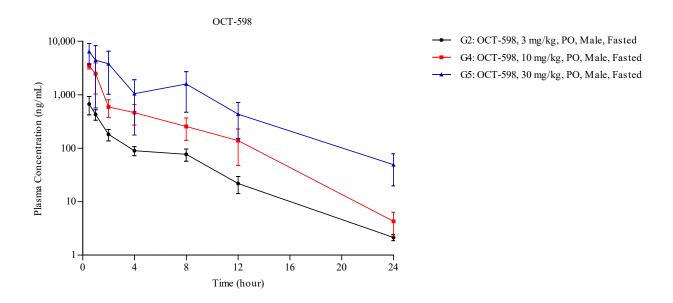


Figure 2. Mean (±SD) Plasma Concentrations-Time Profiles of OCT-598 Following Single Oral Gavage Administration to Fasted Male Rats at 3, 10 and 30 mg/kg

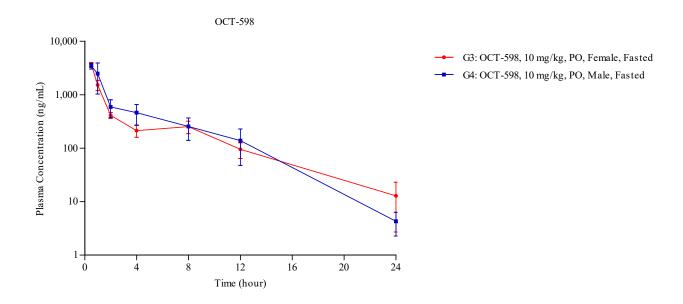


Figure 3. Mean (±SD) Plasma Concentrations-Time Profiles of OCT-598 Following Single Oral Gavage Administration to Fasted Male and Female Rats at 10 mg/kg

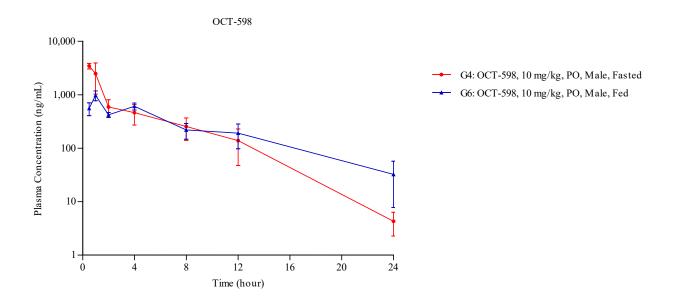


Figure 4. Mean (±SD) Plasma Concentrations-Time Profiles of OCT-598 Following Single Oral Gavage Administration to Fasted and Fed Male Rats at 10 mg/kg

# 10. TABLE

 Table 1. Clinical Observations and Body Weights (Group Summary)

| Group/<br>Treatment | Dose Level (mg/kg) | Dosing<br>Route | Gender | Food<br>Condition | No. of<br>Animals | Clinical Observations | Body Weights (Mean ± SD, g) |
|---------------------|--------------------|-----------------|--------|-------------------|-------------------|-----------------------|-----------------------------|
| G1<br>OCT-598       | 3                  | IV              | Male   | Fed               | 3                 | Appears normal        | $278.0 \pm 12.0$            |
| G2<br>OCT-598       | 3                  | PO              | Male   | Fasted            | 3                 | Appears normal        | $236.2 \pm 10.1$            |
| G3<br>OCT-598       | 10                 | PO              | Female | Fasted            | 3                 | Appears normal        | $173.3 \pm 11.2$            |
| G4<br>OCT-598       | 10                 | РО              | Male   | Fasted            | 3                 | Appears normal        | $239.2 \pm 9.0$             |
| G5<br>OCT-598       | 30                 | РО              | Male   | Fasted            | 3                 | Appears normal        | $230.0 \pm 7.2$             |
| G6<br>OCT-598       | 10                 | РО              | Male   | Fed               | 3                 | Appears normal        | $268.8 \pm 12.4$            |

# 11. APPENDICES

**Appendix 1. Individual Clinical Observations and Body Weights** 

| Group/<br>Treatment | Dose Level (mg/kg) | Dosing<br>Route | Gender    | Food<br>Condition | Animal<br>No. | Clinical Observations | Body Weights (g) |
|---------------------|--------------------|-----------------|-----------|-------------------|---------------|-----------------------|------------------|
|                     |                    |                 |           |                   | M1            | Appears normal        | 258.5            |
|                     |                    |                 |           |                   | M2            | Appears normal        | 294.5            |
| G1                  | 2                  | 13.7            | M.1.      | F. 4              | M3            | Appears normal        | 277.7            |
| OCT-598             | 3                  | IV              | Male      | Fed               | M4            | Appears normal        | 282.3            |
|                     |                    |                 |           |                   | M5            | Appears normal        | 272.4            |
|                     |                    |                 |           |                   | M6            | Appears normal        | 282.7            |
|                     | 3                  | РО              | Male      | Fasted            | M7            | Appears normal        | 234.6            |
|                     |                    |                 |           |                   | M8            | Appears normal        | 247.9            |
| G2                  |                    |                 |           |                   | M9            | Appears normal        | 220.1            |
| OCT-598             |                    |                 |           |                   | M10           | Appears normal        | 231.0            |
|                     |                    |                 |           |                   | M11           | Appears normal        | 245.2            |
|                     |                    |                 |           |                   | M12           | Appears normal        | 238.6            |
|                     |                    |                 |           |                   | F13           | Appears normal        | 179.7            |
|                     |                    | <b>D</b> O      |           | Fasted            | F14           | Appears normal        | 175.1            |
| G3                  | 10                 |                 | PO Female |                   | F15           | Appears normal        | 190.7            |
| OCT-598             | 10                 | PU              |           |                   | F16           | Appears normal        | 159.9            |
|                     |                    |                 |           |                   | F17           | Appears normal        | 170.7            |
|                     |                    |                 |           |                   | F18           | Appears normal        | 163.8            |

|         |    |     |        |         | M19 | Appears normal | 222.2          |
|---------|----|-----|--------|---------|-----|----------------|----------------|
|         |    |     |        |         | M20 | Appears normal | 248.7          |
| G4      | 10 | PO  | M-1-   | Easta d | M21 | Appears normal | 240.7          |
| OCT-598 | 10 | PO  | Male   | Fasted  | M22 | Appears normal | 241.1<br>238.9 |
|         |    |     |        |         | M23 | Appears normal |                |
|         |    |     |        |         | M24 | Appears normal | 243.4          |
|         |    |     |        |         | M25 | Appears normal | 237.8          |
|         |    |     |        |         | M26 | Appears normal | 231.4          |
| G5      | 20 | DO. | N.C. 1 | F 4 1   | M27 | Appears normal | 224.6          |
| OCT-598 | 30 | PO  | Male   | Fasted  | M28 | Appears normal | 227.8          |
|         |    |     |        |         | M29 | Appears normal | 238.2          |
|         |    |     |        |         | M30 | Appears normal | 220.1          |
|         |    |     |        |         | M31 | Appears normal | 254.6          |
|         |    |     |        |         | M32 | Appears normal | 260.1          |
| G6      | 10 | DO. | N.C. 1 | г 1     | M33 | Appears normal | 261.1          |
| OCT-598 | 10 | PO  | Male   | Fed     | M34 | Appears normal | 287.8          |
|         |    |     |        |         | M35 | Appears normal | 275.2          |
|         |    |     |        |         | M36 | Appears normal | 273.9          |

Appendix 2. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 Following Single Intravenous Administration to Fed Male Rats at 3 mg/kg

| Time (1-m) |          | Animal No.     |         | Mann    | CD     | CM (0/) |
|------------|----------|----------------|---------|---------|--------|---------|
| Time (hr)  | Plasma ( | Concentrations | (ng/mL) | Mean    | SD     | CV (%)  |
| 0.002      | M1       | M2             | M3      |         |        |         |
| 0.083      | 3081.74  | 1816.72        | 3034.34 | 2644.27 | 717.07 | 27.1    |
| 0.25       | M4       | M5             | M6      |         |        |         |
| 0.25       | 2000.81  | 2000.48        | 1746.08 | 1915.79 | 146.97 | 7.7     |
| 0.5        | M1       | M2             | M3      |         |        |         |
| 0.5        | 1371.57  | 912.79         | 1266.83 | 1183.73 | 240.41 | 20.3    |
| 1          | M4       | M5             | M6      |         |        |         |
| 1          | 522.73   | 626.75         | 495.88  | 548.45  | 69.13  | 12.6    |
| 2          | M1       | M2             | M3      |         |        |         |
| 2          | 213.46   | 272.08         | 312.31  | 265.95  | 49.71  | 18.7    |
| 4          | M4       | M5             | M6      |         |        |         |
| 4          | 66.79    | 87.36          | 102.71  | 85.62   | 18.02  | 21.1    |
| 0          | M1       | M2             | M3      |         |        |         |
| 8          | 16.92    | 33.77          | 45.87   | 32.18   | 14.54  | 45.2    |
| 2.4        | M4       | M5             | M6      |         |        |         |
| 24         | 3.27     | 1.72           | 2.81    | 2.60    | 0.80   | 30.7    |

Appendix 3. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 Following Single Oral Gavage Administration to Fasted Male Rats at 3 mg/kg

| Time (hm) |          | Animal No.     |         | Mean   | CD     | CV (0/) |
|-----------|----------|----------------|---------|--------|--------|---------|
| Time (hr) | Plasma ( | Concentrations | (ng/mL) | Wiean  | SD     | CV (%)  |
| 0.5       | M7       | M8             | M9      |        |        |         |
| 0.5       | 391.64   | 874.09         | 749.29  | 671.67 | 250.42 | 37.3    |
| 1         | M10      | M11            | M12     |        |        |         |
| 1         | 317.70   | 485.10         | 481.65  | 428.15 | 95.66  | 22.3    |
| 2         | M7       | M8             | M9      |        |        |         |
| 2         | 145.66   | 229.66         | 167.60  | 180.97 | 43.57  | 24.1    |
| 4         | M10      | M11            | M12     |        |        |         |
| 4         | 84.52    | 75.70          | 108.96  | 89.73  | 17.23  | 19.2    |
| 8         | M7       | M8             | M9      |        |        |         |
| δ         | 61.22    | 98.72          | 70.40   | 76.78  | 19.54  | 25.5    |
| 10        | M10      | M11            | M12     |        |        |         |
| 12        | 14.01    | 22.23          | 29.52   | 21.92  | 7.76   | 35.4    |
| 24        | M7       | M8             | M9      |        |        |         |
| 24        | 2.22     | 1.83           | 2.38    | 2.14   | 0.28   | 13.1    |

Appendix 4. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 Following Single Oral Gavage Administration to Fasted Female Rats at 10 mg/kg

| T: (1)    | Animal No. |                |         | Maria   | SD     | CV (0/) |
|-----------|------------|----------------|---------|---------|--------|---------|
| Time (hr) | Plasma (   | Concentrations | (ng/mL) | Mean    | SD     | CV (%)  |
| 0.5       | F13        | F14            | F15     |         |        |         |
| 0.5       | 3507.42    | 3546.25        | 4051.96 | 3701.88 | 303.80 | 8.2     |
| 1         | F16        | F17            | F18     |         |        |         |
| 1         | 1149.70    | 1667.08        | 1767.72 | 1528.17 | 331.60 | 21.7    |
| 2         | F13        | F14            | F15     |         |        |         |
| 2         | 429.58     | 353.40         | 454.72  | 412.56  | 52.76  | 12.8    |
| 4         | F16        | F17            | F18     |         |        |         |
| 4         | 175.01     | 190.82         | 272.43  | 212.75  | 52.29  | 24.6    |
| 0         | F13        | F14            | F15     |         |        |         |
| 8         | 193.48     | 241.21         | 325.09  | 253.26  | 66.63  | 26.3    |
| 10        | F16        | F17            | F18     |         |        |         |
| 12        | 63.25      | 97.43          | 125.15  | 95.28   | 31.00  | 32.5    |
| 24        | F13        | F14            | F15     |         |        |         |
| 24        | 3.04       | 12.23          | 23.33   | 12.87   | 10.16  | 79.0    |

BLQ < LLOQ (1 ng/mL)

Appendix 5. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 Following Single Oral Gavage Administration to Fasted Male Rats at 10 mg/kg

| Time (hm) |          | Animal No.     |         | Mean    | SD      | CV (0/) |
|-----------|----------|----------------|---------|---------|---------|---------|
| Time (hr) | Plasma ( | Concentrations | (ng/mL) | Mean    | SD      | CV (%)  |
| 0.5       | M19      | M20            | M21     |         |         |         |
| 0.5       | 3410.52  | 3880.34        | 3042.86 | 3444.57 | 419.78  | 12.2    |
| 1         | M22      | M23            | M24     |         |         |         |
| 1         | 2149.51  | 1223.88        | 4075.58 | 2482.99 | 1454.80 | 58.6    |
| 2         | M19      | M20            | M21     |         |         |         |
| 2         | 697.94   | 728.99         | 345.54  | 590.82  | 212.99  | 36.0    |
| 4         | M22      | M23            | M24     |         |         |         |
| 4         | 301.41   | 411.97         | 674.93  | 462.77  | 191.87  | 41.5    |
| 8         | M19      | M20            | M21     |         |         |         |
| 8         | 179.45   | 386.04         | 197.90  | 254.46  | 114.32  | 44.9    |
| 10        | M22      | M23            | M24     |         |         |         |
| 12        | 165.68   | 37.19          | 212.68  | 138.52  | 90.85   | 65.6    |
| 24        | M19      | M20            | M21     |         |         |         |
| 24        | 5.62     | 5.33           | 1.97    | 4.30    | 2.03    | 47.2    |

Appendix 6. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 Following Single Oral Gavage Administration to Fasted Male Rats at 30 mg/kg

| Time (hr) |          | Animal No.     |         | Mean    | SD      | CV (0/) |  |
|-----------|----------|----------------|---------|---------|---------|---------|--|
| Time (hr) | Plasma ( | Concentrations | (ng/mL) | Mean    | SD      | CV (%)  |  |
| 0.5       | M25      | M26            | M27     |         |         |         |  |
| 0.5       | 4103.79  | 9343.71        | 5786.04 | 6411.18 | 2675.31 | 41.7    |  |
| 1         | M28      | M29            | M30     |         |         |         |  |
| 1         | 1084.43  | 8647.34        | 3544.14 | 4425.30 | 3857.69 | 87.2    |  |
| 2         | M25      | M26            | M27     |         |         |         |  |
| 2         | 1443.66  | 6831.94        | 3082.55 | 3786.05 | 2762.17 | 73.0    |  |
| 4         | M28      | M29            | M30     |         |         |         |  |
| 4         | 367.78   | 2024.90        | 742.56  | 1045.08 | 868.99  | 83.2    |  |
| 0         | M25      | M26            | M27     |         |         |         |  |
| 8         | 599.71   | 2799.60        | 1367.10 | 1588.80 | 1116.58 | 70.3    |  |
| 12        | M28      | M29            | M30     |         |         |         |  |
| 12        | 119.47   | 516.81         | 663.83  | 433.37  | 281.61  | 65.0    |  |
| 2.4       | M25      | M26            | M27     |         |         |         |  |
| 24        | 29.72    | 83.09          | 34.90   | 49.24   | 29.43   | 59.8    |  |

Appendix 7. Individual Plasma Concentrations and Pharmacokinetic Parameters of OCT-598 Following Single Oral Gavage Administration to Fed Male Rats at 10 mg/kg

| Time a (law) |          | Animal No.     |         | Maan   | CD     | CV (0/) |
|--------------|----------|----------------|---------|--------|--------|---------|
| Time (hr)    | Plasma ( | Concentrations | (ng/mL) | Mean   | SD     | CV (%)  |
| 0.5          | M31      | M32            | M33     |        |        |         |
| 0.3          | 629.02   | 384.17         | 657.21  | 556.80 | 150.16 | 27.0    |
| 1            | M34      | M35            | M36     |        |        |         |
| 1            | 884.90   | 833.95         | 1209.26 | 976.03 | 203.58 | 20.9    |
| 2            | M31      | M32            | M33     |        |        |         |
| 2            | 372.01   | 439.69         | 451.89  | 421.19 | 43.03  | 10.2    |
| 4            | M34      | M35            | M36     |        |        |         |
| 4            | 506.91   | 641.35         | 677.86  | 608.71 | 90.03  | 14.8    |
| 0            | M31      | M32            | M33     |        |        |         |
| 8            | 254.27   | 267.41         | 136.85  | 219.51 | 71.89  | 32.7    |
| 10           | M34      | M35            | M36     |        |        |         |
| 12           | 261.97   | 85.15          | 226.20  | 191.11 | 93.49  | 48.9    |
| 2.4          | M31      | M32            | M33     |        |        |         |
| 24           | 60.87    | 19.45          | 16.91   | 32.41  | 24.68  | 76.2    |


# Appendix 8. Certificate of Analysis (COA)

| contr  | EM        | edl  | Certific                                | ate              | of Analys  |  | facturer<br>iss and<br>Contact       |  | QF0416<br>Asymci<br>No. 71, 7th Avenue,  | nem Life Science (<br>TEDA Tianjin, 300 | 09/2021<br>Tianjin) Co., Lt<br>1457, P.R. Chir<br>(0)22 6625288 |
|--|-----------|--|---|------------------|--|--|--------------------------------------|--|--|---|---|
|  |           |  |   |                  |  | C                                      | DA No.                               |  |  | TJ2-CP-20                               | 2305-074.0  |
| Basic information  |           |  |   |                  |  |  |                                      |  |  |   |   |
| Asymchem   |           | 404704   | Client product                          |                  | 007.500  | CAS No.                                |                                      |  |  |   | _   |
| Asymchem   | CCL No.   | January Company  |   |                  |  |  | 51755/2                              |  | . 0  |   | 11  |
| Chemic   | al Name   | (2R,4r,6R)-6-(4  | 4-((3'-fluoro-5'-me<br>nene-3-carboxami | ethoxy<br>ido)sp | -[1,1'-biphenyl]-4-y<br>iro[3,3]heptane-2-c  | .)methyl)-2<br>arboxylic a             | ,5-<br>icid                          |  |  | X XIVIX                                 | OH  |
| Molecular  | Formula   | C29H30FNO4S  |   |                  | Molecular Weight   |  |                                      |  | S  | . ~ ~                                   | `H  |
|  |           |  |   |                  |  |  |                                      |  |  |   |   |
|  | Lot No.   | CPo131731-0  | 1-08-01                                 |                  | Batch size   | 1.740 kg                               |                                      |  | 1  |   |   |
| Date of mar  | ufacture  | May 27, 2023   |   |                  | Testing date   | May 30, 2                              | 023                                  |  | 5  |   |   |
| Release testing  |           | Asymchem Life Sc   | cience (Tianjin) Co., Lt                | d.               |  |  |                                      |  |  |   |   |
| -  |           |  | e, TEDA Tianjin, 3004                   | 57, P.R.         | China  |  |                                      |  |  |   |   |
| Storage  | condition | 5°C (2°C ~8°C  | ;)                                      |                  |  |  |                                      |  | E  |   |   |
| Re   | test date | May 26, 2024   |   |                  |  |  |                                      |  | 0-   | -                                       |   |
| Testing results  |           |  |   |                  |  |  |                                      |  |  |   |   |
|  |           |  |   | _                |  |  |                                      |  | NOTEBOOK   |   |   |
| TEST ITEM  |           | Target Quali   | ity*                                    | A                | NALYTICAL RESU   | LTS                                    | TEST                                 | METHODS  | REFERENCE  | Pass / Fail                             | Commer  |
| Appearance   |           | Report   |   |                  | Off-white solid  |  | AM-                                  | -001.03  | TJ2-N01173   | -                                       | -   |
| Appearance-  |           | December   |   | 2000 00 00 00    |  |  |                                      |  | page 256   |   |   |
| solution   |           | Report   |   |                  |  | 100                                    |                                      |  | T.I2-N03114  |   |   |
|  |           | Papart   |   | L                | ight yellow clear lic  | quid                                   | AM-131                               | 731-018.01   | TJ2-N03114<br>page 20-21   | -                                       | \ =   |
| Identification<br>(FTIR)   |           | Report   |   | L                | ight yellow clear lic<br>Refer to attachme   | ******                                 |                                      | 731-018.01   | page 20-21<br>TJ2-N00876   | -                                       | -   |
| (FTIR)<br>Identification   |           |  | ucture                                  |                  | Refer to attachme  | nt                                     | AM-                                  | -023.04  | page 20-21   | -                                       |   |
| (FTIR)<br>Identification<br>(1H-NMR)   | A         | Conforms to str  |   |                  | Refer to attachme  | nt<br>ure                              | AM-                                  |  | page 20-21<br>TJ2-N00876<br>page 201-202<br>TJ2-N01120<br>page 53-54   |   | -   |
| (FTIR)<br>Identification   | A         |  |   |                  | Refer to attachme  | nt<br>ure                              | AM-                                  | -023.04  | page 20-21<br>TJ2-N00876<br>page 201-202<br>TJ2-N01120<br>page 53-54<br>TJ2-N02137   | -                                       |   |
| (FTIR) Identification (1H-NMR) Identification (HPLC) Identification  | RT ma     | Conforms to stri<br>atches to that o<br>standard<br>ctrum (diode an  | of reference                            | RT m             | Refer to attachme Conforms to structivatches to that of restandard UV spectrum match   | nt<br>ure<br>iference                  | AM-<br>AM-131                        | 023.04<br>014.04<br>731-016.02                                       | page 20-21<br>TJ2-N00876<br>page 201-202<br>TJ2-N01120<br>page 53-54<br>TJ2-N02137<br>page 95-102<br>TJ2-N01903  | -<br>Pass<br>Pass                       | -   |
| (FTIR) Identification (1H-NMR) Identification (HPLC) Identification (UV)                                     | RT ma     | Conforms to structure (diode arms of reference   | ray) matches                            | RT m             | Refer to attachme Conforms to struct natches to that of re standard UV spectrum match hat of reference sta   | nt<br>ure<br>sference<br>nes<br>andard | AM-<br>AM-131<br>AM-131              | 023.04<br>014.04<br>731-016.02<br>731-024.01                         | page 20-21 TJ2-N00876 page 201-202 TJ2-N01120 page 53-54 TJ2-N02137 page 95-102 TJ2-N01903 page 236-240  | -<br>Pass                               | -   |
| (FTIR) Identification (1H-NMR) Identification (HPLC) Identification  | RT ma     | Conforms to stri<br>atches to that o<br>standard<br>ctrum (diode an  | ray) matches                            | RT m             | Refer to attachme Conforms to structivatches to that of restandard UV spectrum match   | nt ure eference nes indard rence       | AM-<br>AM-131<br>AM-131              | 023.04<br>014.04<br>731-016.02                                       | page 20-21<br>TJ2-N00876<br>page 201-202<br>TJ2-N01120<br>page 53-54<br>TJ2-N02137<br>page 95-102<br>TJ2-N01903<br>page 236-240<br>TJ2-N01297<br>page 193-195  | -<br>Pass<br>Pass                       | -   |
| (FTIR) Identification (1H-NMR) Identification (HPLC) Identification (UV) Polymorph Assay (% w/w)             | RT ma     | Conforms to structure atches to that of standard ctrum (diode an eat of reference institute to the reference in the r | ray) matches                            | RT m             | Refer to attachme Conforms to structivatches to that of restandard UV spectrum matchat of reference standards to the reference st | nt ure eference nes indard rence       | AM-<br>AM-131<br>AM-131<br>AM-       | 023.04<br>014.04<br>731-016.02<br>731-024.01                         | page 20-21<br>TJ2-N00876<br>page 201-202<br>TJ2-N01120<br>page 53-54<br>TJ2-N02137<br>page 95-102<br>TJ2-N01903<br>page 236-240<br>TJ2-N01297<br>page 193-195<br>TJ2-N02137                              | Pass Pass Pass                          | -   |
| (FTIR) Identification (1H-NMR) Identification (HPLC) Identification (UV) Polymorph Assay (% w/w) As is basis | RT ma     | Conforms to stra<br>atches to that o<br>standard<br>ctrum (diode an<br>nat of reference<br>ns to the refere<br>(Form A)  | ray) matches                            | RT m             | Refer to attachme Conforms to structivatches to that of restandard UV spectrum matches of reference standard (Form A 96.6%)  | nt ure eference nes indard rence       | AM-<br>AM-131<br>AM-131<br>AM-131    | 023.04<br>014.04<br>731-016.02<br>731-024.01<br>029.08<br>731-016.02 | page 20-21<br>TJ2-N00876<br>page 201-202<br>TJ2-N01120<br>page 53-54<br>TJ2-N02137<br>page 95-102<br>TJ2-N01903<br>page 236-240<br>TJ2-N01297<br>page 193-195  | Pass Pass Pass Pass                     | -   |
| (FTIR) Identification (1H-NMR) Identification (HPLC) Identification (UV) Polymorph Assay (% w/w)             | RT ma     | Conforms to stratches to that of standard ctrum (diode an nat of reference ins to the reference (Form A)   | ray) matches                            | RT m             | Refer to attachment Conforms to structionatches to that of restandard UV spectrum materials of reference standard (Form A sta | nt ure eference nes indard rence       | AM-<br>AM-131<br>AM-131<br>AM-131    | 023.04<br>014.04<br>731-016.02<br>731-024.01<br>029.08               | page 20-21<br>TJ2-N00876<br>page 201-202<br>TJ2-N01120<br>page 53-54<br>TJ2-N02137<br>page 95-102<br>TJ2-N01903<br>page 236-240<br>TJ2-N01297<br>page 193-195<br>TJ2-N02137<br>page 95-102               | Pass Pass Pass                          | -   |
| (FTIR) Identification (1H-NMR) Identification (HPLC) Identification (UV) Polymorph Assay (% w/w) As is basis | RT ma     | Conforms to stra<br>atches to that o<br>standard<br>ctrum (diode an<br>nat of reference<br>ns to the refere<br>(Form A)  | ray) matches                            | RT m             | Refer to attachme Conforms to structivatches to that of restandard UV spectrum matches of reference standard (Form A 96.6%)  | nt ure eference nes indard rence       | AM-131<br>AM-131<br>AM-131<br>AM-131 | 023.04<br>014.04<br>731-016.02<br>731-024.01<br>029.08<br>731-016.02 | page 20-21<br>TJ2-N00876<br>page 201-202<br>TJ2-N01120<br>page 53-54<br>TJ2-N02137<br>page 95-102<br>TJ2-N01903<br>page 236-240<br>TJ2-N01297<br>page 193-195<br>TJ2-N02137<br>page 95-102<br>TJ2-N00867 | Pass Pass Pass Pass                     | -   |

Asymchem Confidential Information

| ASYMCHE   | Med   | Certific          | cate of Analysis   | Manufactu<br>Address ar<br>Cont | nd              | QF04168.V06 (03/09/2021<br>Asymchem Life Science (Tianjin) Co.,<br>No. 71, 7th Avenue, TEDA Tianjin, 300457, P.R. Cl<br>Tel: +86 (0)22 66252. |                      |        |  |  |
|---|---|-------------------|--|---------------------------------|-----------------|---|----------------------|--------|--|--|
| 0 011616  |   |                   |  |                                 | lo.             | - 1 2 2   | TJ2-CP-202305-074.02 |        |  |  |
| TEST ITEM   | Targe   | t Quality*        | ANALYTICAL RESULTS   | 3                               | TEST METHODS    | NOTEBOOK<br>REFERENCE   | Pass / Fail          | Commen |  |  |
| Any single<br>unspecified<br>impurity<br>(% area) | Report RRT and % area for all<br>impurities at or above 0.05%<br>Target 0.5% for maximum for a<br>single impurity |                   | RRT 0.42: 0.11%<br>RRT 0.93: 0.12%<br>RRT 1.22: 0.97%<br>RRT 1.33: 0.15%<br>RRT 1.37: 0.44%<br>RRT 1.38: 0.19% | A                               | M-131731-016.02 | TJ2-N02137<br>page 95-102   | -                    | -      |  |  |
| Total impurity (% area)                           |   | eport<br>about 3% | 2.7%   |                                 |                 |   | -                    | -      |  |  |
| Chiral purity (% area)                            | NLT 99.0%   |                   | 99.6%  | A                               | M-131731-021.01 | TJ2-N01324<br>page 117-122  | Pass                 | -      |  |  |
| Residual solvent                                  |   |                   |  |                                 |                 |   |                      |        |  |  |
| Dichloromethane                                   | NMT   | 600ppm            | <60 ppm  |                                 |                 |   | Pass                 | -      |  |  |
| 2-Methyl THF                                      | NMT :   | 5000ppm           | ND (not detected)  |                                 |                 |   | Pass                 | -      |  |  |
| n-Heptane   | NMT :   | 5000ppm           | <500 ppm   |                                 |                 |   | Pass                 | -      |  |  |
| Methanol  | NMT :   | 3000ppm           | ND (not detected)  | ND (not detected)               |                 |   | Pass                 | 2:     |  |  |
| THF   | NMT   | 720ppm            | ND (not detected)  |                                 | M-131731-025.01 | TJ2-N01300<br>page 238-249  | Pass                 | 2      |  |  |
| MTBE  | NMT :   | 5000ppm           | ND (not detected)  | A                               | W-131/31-025.01 |   | Pass                 | -      |  |  |
| Ethyl acetate                                     | NMT :   | 5000ppm           | 735 ppm  |                                 |                 |   | Pass                 | -      |  |  |
| NMP   | NMT   | 530ppm            | ND (not detected)  |                                 |                 |   | Pass                 | -      |  |  |
| Triethylamine                                     | R   | eport             | <500 ppm   |                                 |                 |   | -                    | 35     |  |  |
| DIPEA   | R   | eport             | <150 ppm   |                                 |                 |   | -                    | -      |  |  |
| Water Content<br>(% w/w)                          | R   | eport             | 0.14%  |                                 | AM-003.08       | TJ2-N01173<br>page 257-259  | -                    | -      |  |  |
| Residue on<br>ignition                            | NM  | Γ 0.5%            | 0.0% (0.01%)   |                                 | USP <281>       | TJ2-N02354<br>page 77-84  | Pass                 | -      |  |  |
|   | Cd  | ≤0.5 ppm          | ND, <lod, lod="5" ppb<="" td=""><td></td><td></td><td>T.J2-2023-</td><td>Pass</td><td>-</td></lod,>            |                                 |                 | T.J2-2023-  | Pass                 | -      |  |  |
| Elemental impurities (ppm)                        | Pb  | ≤0.5 ppm          | 0.0 ppm (0.005 pm)   |                                 | AM-066.05       | 00309   | Pass                 | 1-     |  |  |
|   | As ≤1.5 ppm   |                   | ND, <lod, lod="30" ppb<="" td=""><td colspan="2"></td><td>page 1-10</td><td>Pass</td><td>-</td></lod,>         |                                 |                 | page 1-10   | Pass                 | -      |  |  |



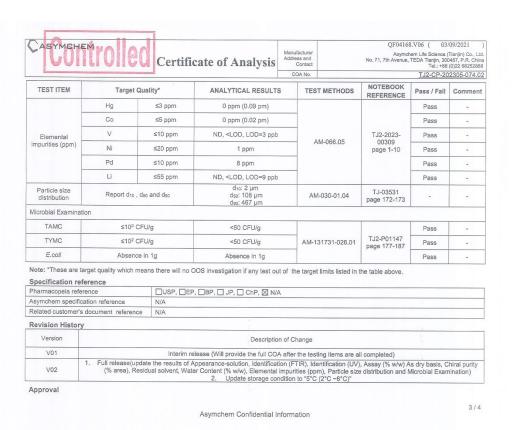

QF04168.V06 ( 03/09/2021 ASYMCHEM | Asymchem Life Science (Tianjin) Co., Ltd.
No. 71, 7th Avenue, TEDA Tianjin, 300457, P.R. China
Tel.: +86 (0)22 66252888 Certificate of Analysis TJ2-CP-202305-074.02 MA Signature/Date: 6106 60 170 Reviewed by Lab Director Ma 57/03/20rs Approved by QA QA Conclusion Pass QA Hold NA Signature/Date: 07/03/2016

4/4

Asymchem Confidential Information

iei

# Appendix 9. Analytical Method

# 1. METHOD SUMMARY

| Analyte / Metabolite                        | OCT-598 (Lot no. CPo131731-01-08-01)   |
|---|--|
|   | `  |
| Internal Standard                           | OCT-598-D3 (Lot no. P-0003800-002)   |
| Species (Strain) / Matrix                   | Rat, Sprague-Dawley / Plasma   |
| Anticoagulant                               | Sodium heparin   |
| Sample Volume                               | 20 μL  |
| Calibration Range                           | 1 to 1000 ng/mL  |
| Regression Type / Weighting                 | Linear / 1/x <sup>2</sup>  |
| Extraction                                  | Acetonitrile protein precipitation   |
| LC-MS/MS                                    | 1290 Infinity II, Agilent / Triple Quad 4500, Sciex (QuBEST Instrument ID: 4500-1) |
| Sample Processing Temperature               | Room temperature   |
| Sample Storage Temperature                  | -80°C  |
| Special Storage / Treatment<br>Requirements | None   |

# 2. LC CONDITION

| Column                    | ZORBAX RRHD Eclipse Plus C8 (2.1 x 100 mm, 1.8 μm), Agilent (ID: C8-119) |                     |       |       |  |  |  |
|---------------------------|--|---------------------|-------|-------|--|--|--|
| Column Temperature        | 40°C   | 40°C                |       |       |  |  |  |
| Autosampler Temperature   | 15°C   |                     |       |       |  |  |  |
| Injection Volume          | 2 μL   |                     |       |       |  |  |  |
| Needle Wash Solution      | Isopropyl alcohol:Ultrapure water (UPW), 9:1 (v/v) (Flush Port, 30sec.)  |                     |       |       |  |  |  |
| Mobile Phase A            | 10 mM Ammor  | nium acetate in UPV | V     |       |  |  |  |
| Mobile Phase B            | 0.1% Formic ac   | cid in acetonitrile |       |       |  |  |  |
| Mobile Phase Time Program | Time (min)   | Flow (mL/min)       | A (%) | B (%) |  |  |  |
|                           | 0  | 0.4                 | 35    | 65    |  |  |  |
|                           | 2.5  | 0.4                 | 35    | 65    |  |  |  |

# 3. MASS SPECTROMETER CONDITION

| Parameter             | OCT-598       | OCT-598-D3 (IS) |  |  |  |  |
|-----------------------|---------------|-----------------|--|--|--|--|
| MRM (m/z)             | 508.2 > 306.0 | 511.2 > 306.0   |  |  |  |  |
| DP (V)                | 121           | 126             |  |  |  |  |
| EP (V)                | 10            | 10              |  |  |  |  |
| CE (V)                | 31            | 31              |  |  |  |  |
| CXP (V)               | 10            | 12              |  |  |  |  |
| Polarity              | ESI Positive  |                 |  |  |  |  |
| CUR (psi)             | 30            |                 |  |  |  |  |
| GS1 (psi)             |               | 30              |  |  |  |  |
| GS2 (psi)             |               | 40              |  |  |  |  |
| CAD (psi)             | F             | High            |  |  |  |  |
| Ion spray voltage (V) | 5             | 5500            |  |  |  |  |
| Ion spray temp. (°C)  | :             | 550             |  |  |  |  |
| Nebulizing gas        | Nit           | trogen          |  |  |  |  |
| Data processing       | Anal          | yst 1.7.0       |  |  |  |  |

# 4. EXTRACTION PROCEDURE

| No. | Procedure  |
|-----|--|
| 1   | The samples exceeded the ULOQ concentration of the calibration curve at each         |
|     | concentration were used after 10, 20, 30 fold dilution by adding blank sample.       |
| 2   | 20 μL of study sample was transferred into a microtube. However, tubes containing    |
|     | 20 μL of standards (CAL or QC sample) were prepared.                                 |
| 3   | 300 μL of IS working solution (ISW3, 5 ng/mL) was added into the microtube.          |
|     | However, for DB sample, added 300 μL of acetonitrile instead of IS working solution. |
| 4   | The mixture was mixed at 2000 rpm for 2 minutes using vortex mixer, and then         |
| 4   | centrifuged at 15000 RCF at 4°C for 5 minutes.                                       |
| 5   | 200 μL of supernatant was transferred into a sample vial.                            |